

sealed. sterile. secure.

For complete cleanroom & isolator confidence

Thermo Scientific Triple Wrapped Irradiated Plates

Industry-leading performance inside and out

Designed for your most regulated environments, Thermo Scientific™ Triple Wrapped Irradiated Plates feature the latest in quality-assurance technology. From advanced layered protection to enhanced performance, it's the superior choice for environmental monitoring.

Confidence in results

- VHP exposure indicator confirms the integrity of the wrapping material and seals during VHP exposure, ensuring optimum media fertility

Contamination risk reduction


- Protective layers enhance moisture control by incorporating novel elements that permit the passage and absorption of excess moisture, reducing the risk of contamination

Cost-effective solutions

- Improved shelf life and room temperature storage minimize new batch validation and refrigeration costs



There is no room to compromise when it comes to your regulatory compliance.

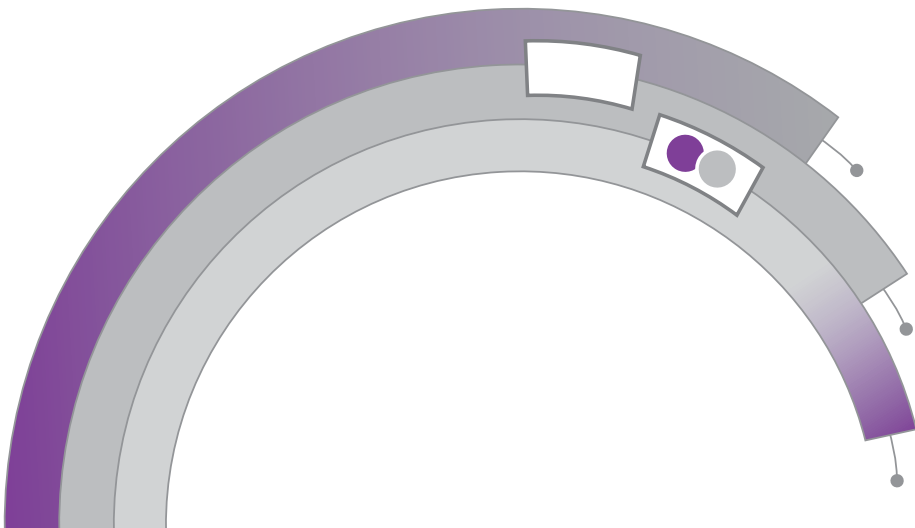


Thermo
SCIENTIFIC

Committed to safety

A range of high quality, Thermo Scientific™ prepared culture media products, which are ideal for environmental monitoring in the pharmaceutical industry. These include 90mm deep fill plates for passive and active air sampling applications. Containing Tryptone Soya Agar, a general purpose medium, or Sabouraud Dextrose Agar, ideal for the growth of yeasts and moulds, they are available with or without neutralisers. The plates are irradiated and triple-wrapped in transparent packaging. Includes an exposure indicator to reduce contamination risk in Vapour-phase Hydrogen Peroxide (VHP) sterilisation procedures.

55mm contact plates are designed for the microbiological testing of surfaces and personnel, and the new format makes them even easier to store, handle and use. Grip-lock lid helps to prevent accidental exposure of the agar and ensure the integrity of the medium during handling and transport. Automated manufacturing procedures ensure a smooth, evenly distributed layer of agar on each plate. A grid on the base of the plate subdivides the growth area to aid accurate colony counts.

All products undergo quantitative quality control testing, using a low-dose inoculum to ensure optimum performance. The range provides high quality culture media for total aerobic microbial counts, coliform bacteria counts and yeast/mould counts, for use in designated clean areas.

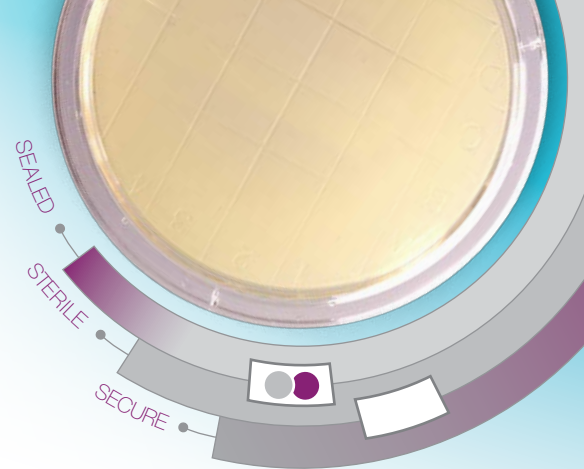

Sealed

- Transparent easy open wrapping materials
- Primary wrap PET/PE layer with breathable patch
- Secondary wrap PET/SiOx/PE layer impervious to VHP with immobilised moisture control sachet
- Tertiary wrap PET/PE layer dust cover
- 10mm diameter hanging hole in end seal of primary and secondary wraps

Sterile

- Sterility assurance level SAL 10⁻⁵
- Integrated irradiation color indicator

Secure

- Media storage temperature range 2–25° C
- Enhanced product shelf life (see product specification)
- Integrated VHP exposure indicator
- Batch specific certificates of origin for animal components
- Quantitative quality control testing to ensure optimum performance in accordance with BP/EP/JP/USP

Thermo Scientific Triple Wrapped Irradiated Plates

Order Information		
Description	Packaging (2x5)	Global Product Code
Triple Wrap Sterile Pack w/VHP Indicator Tryptone Soya Agar	Settling Plates 90mm	P05500B
Triple Wrap Sterile Pack w/VHP Indicator Tryptone Soya Agar w/Lecithin, Polysorbate 80, Sodium Thiosulphate, L-Histidine	Settling Plates 90mm	P05501B
Triple Wrap Sterile Pack w/VHP Indicator Sabouraud Dextrose Agar	Settling Plates 90mm	P05502B
Triple Wrap Sterile Pack w/VHP Indicator Sabouraud Dextrose Agar w/Lecithin, Polysorbate 80, Sodium Thiosulphate, L-Histidine	Settling Plates 90mm	P05503B
Triple Wrap Sterile Pack w/VHP Indicator Tryptone Soya Agar	Contact Plates 55mm	P05510D
Triple Wrap Sterile Pack w/VHP Indicator Tryptone Soya Agar w/Lecithin, Polysorbate 80, Sodium Thiosulphate, L-Histidine	Contact Plates 55mm	P05511D
Triple Wrap Sterile Pack w/VHP Indicator Sabouraud Dextrose Agar	Contact Plates 55mm	P05512D
Triple Wrap Sterile Pack w/VHP Indicator Sabouraud Dextrose Agar w/Lecithin, Polysorbate 80, Sodium Thiosulphate, L-Histidine	Contact Plates 55mm	P05513D

thermoscientific.com/triplewrap

© 2013 Thermo Fisher Scientific. All rights reserved. All trademarks are property of Thermo Fisher Scientific and its subsidiaries.

Contact Information:

International
+44 (0) 1256 841144
oxoid.info@thermofisher.com

USA
+1 800 255 6730
microbiology@thermofisher.com

LT2098A Rev2
October 2013

Thermo
SCIENTIFIC

Part of Thermo Fisher Scientific